

HD80 HAND CONTROL

Features:

- Colour: grey (RAL 7035)
- Protection class: IPx4
- Service indicator that blinks:
 - Every 12 months (norm EN10535 recommendation)/ after 8000 cycles (whichever comes first)
(Via a special Hand control it is possible to reprogramme the CBJ-Care and set other service intervals)
- Battery indicators:
 - Green; battery fully charged (approx. 100-50% capacity remaining)
 - Green; battery about half full (approx. 50-25% capacity remaining)
 - Yellow; Battery needs charging (less than 25% capacity remaining. Incl. audio signal when the hand control is activated)
- Overload indicator:
 - Will blink when the max learned weight is exceeded (learn mode for current cut off value)
- Only compatible with CBJ-C.

Options:

- Customised front cover
- Available with hook or magnet

Usage:

- Easy resetting of service indicator via HD80



MEDLINE[®]
IMPROVING EFFICIENCY
CARELINE[®]
IMPROVING EFFICIENCY



The HD80 is a new handset with an optimized ergonomic design and functions that are activated via "dome" buttons.

The front cover is especially designed for patient lifters and JUMBO Care (CBJ-C) with service monitor.

The HD80 has a logical front cover drawing, a 3 stage battery indication and a service indicator on the hand control. Furthermore the hand control has an overload indicator blinking if the max intended load for the lift is reached.

The informing diodes ensure relevant information always at hand.

**Dimensions (mm):
HB81 - HB85**

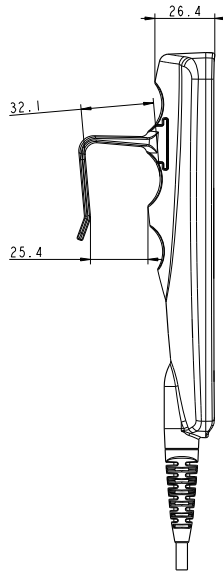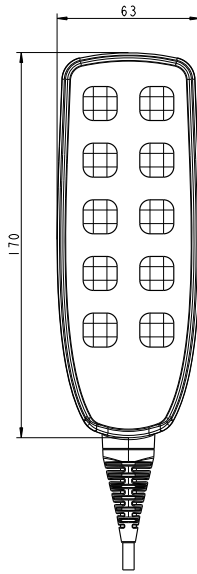


With hook



With magnet

Front covers and ordering numbers:

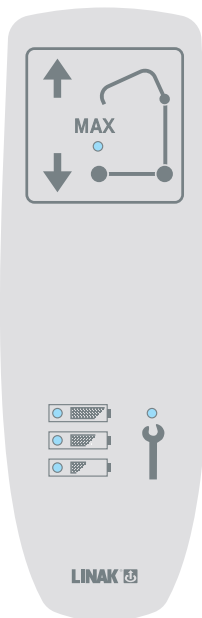

Ordering number:
HD8x008-00 (With hook)
HD8x010-00 (With magnet)



Ordering number:
HD8x009-00 (with hook)
HD8X011-00 (with magnet)

Terms of use

The user is responsible for determining the suitability of LINAK products for specific application. LINAK takes great care in providing accurate and up-to-date information on its products. However, due to continuous development in order to improve its products, LINAK products are subject to frequent modifications and changes without prior notice. Therefore, LINAK cannot guarantee the correct and actual status of said information on its products. While LINAK uses its best efforts to fulfil orders, LINAK cannot, for the same reasons as mentioned above, guarantee the availability of any particular product. Therefore, LINAK reserves the right to discontinue the sale of any product displayed on its website or listed in its catalogues or other written material drawn up by LINAK. All sales are subject to the Standard Terms of Sale and Delivery for LINAK. For a copy hereof, please contact LINAK.